



# OMNIPLUS OSM200

**DENTISTRY**

**SURGICAL & CLINICAL  
DOCUMENT IN 4K**



**NATIONWIDE COVERAGE  
SEE BETTER**

**PRESCOTT'S, INC.**  
*Surgical Microscope Specialists*

# OMNIPLUS OSM200

*For better results, every day.*

Incorporate everything needed to fully concentrate on the procedure at hand.

## TECHNICAL DATA

<b>Magnification</b>	5-step apochromatic magnification changer
<b>Focus</b>	Objective lens 200-350mm options Fine focus or varioscope
<b>Main Lens</b>	f=200-350mm options Apochromatic
<b>Binocular</b>	0-180 inclinable with ergonomic angled wedge; 10x or 12.5x eyepieces
<b>Illumination</b>	60w LED (5500K) paired with liquid light guide; optimal transmission; orange filter
<b>Video</b>	4K documentation, still and live image
<b>Suspension</b>	Rolling stand, wall and ceiling mount
<b>Asepsis</b>	Drapes, drape ring adapter, sterilizable caps
<b>Accessories</b>	Extension arm, rotational plate, MORA, onboard monitor

## SAVE ON SPACE

With rolling stand, wall mount, and ceiling mount options.



## DENTAL APPLICATIONS

- Endodontics
- Restorative
- Aesthetic
- Periodontics
- Implantology
- Educational use

## OPTIMAL CLARITY SUPERIOR FOCUS & MAGNIFICATION

- 0-180 inclinable binoculars
- Ergonomic wedge
- 10x and 12.5x eyepieces
- High precision optics with broadband anti-reflexive coating
- Apochromatic 5-step magnification changer
- Lens options: 250-350mm focusable and varioscope

## SEE BETTER LED ILLUMINATION & LLG

- Superior light output allows for better viewing of deep and narrow canals
- Liquid light guide (LLG) guarantees maximum light transmission
- 60w LED ideal for UV curing
- Color Rendering Index rating of 92
- Resembles natural daylight at 5700K color temperature
- 50,000 hours lifespan
- Easily adjust light intensity
- Accepts multiple cable sizes

## DOCUMENT IN 4K STILL & LIVE IMAGE

- Integrated 4k camera
- Education and lecture uses
- SD memory card
- Wireless remote option


## CUSTOMIZABLE TO MEET YOUR NEEDS

- Wall mount, ceiling mount, rolling stand options
- Long or short swing arm
- Integrated or non-integrated video
- Onboard monitor with mount
- Extension arm
- MORA attachments
- F170/F260 binoculars



NATIONWIDE COVERAGE

PRESCOTT'S, INC.  
*Surgical Microscope Specialists*

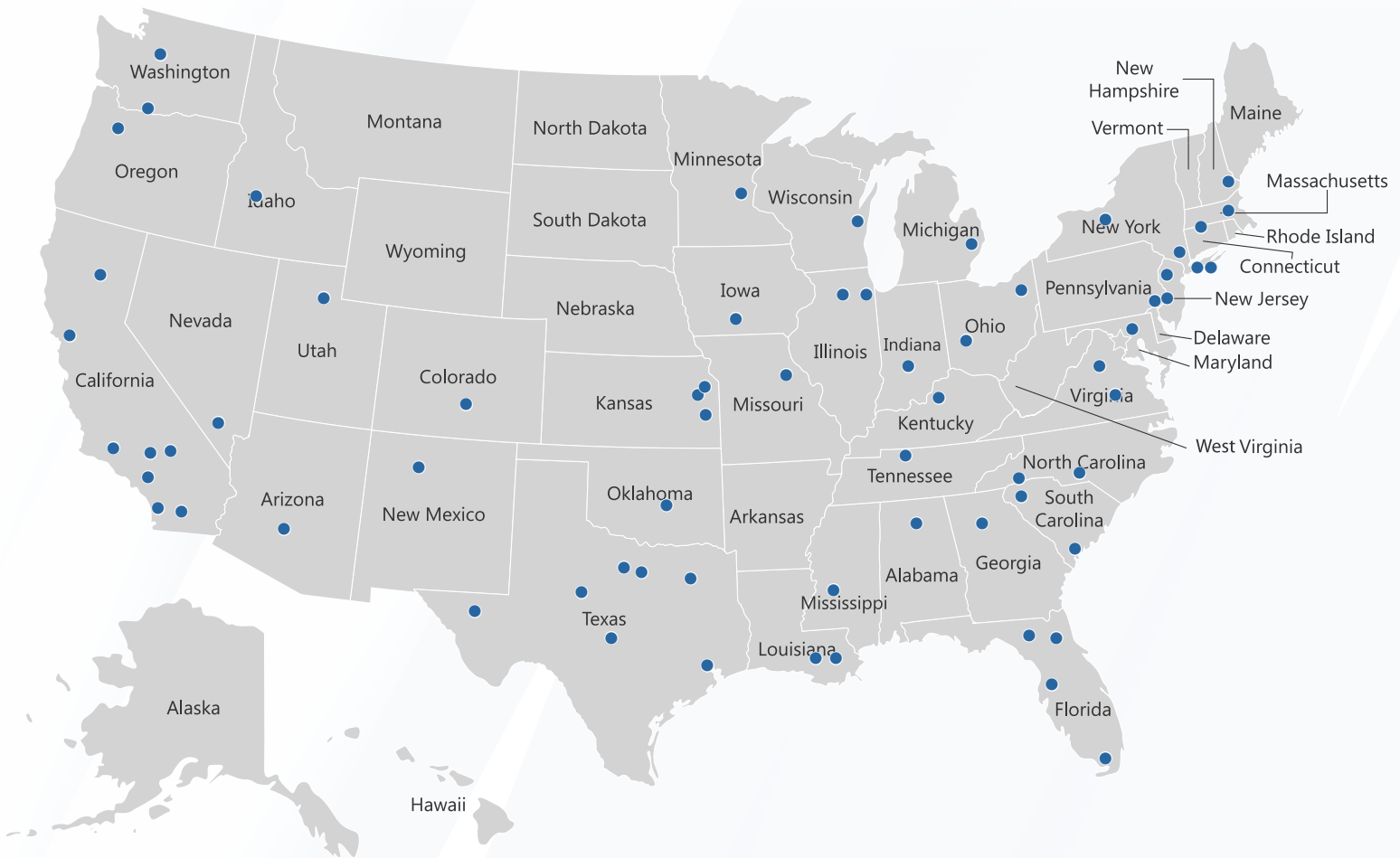

## CONTACT PRESCOTT'S TODAY TO GET IN TOUCH WITH YOUR EXCLUSIVE REPRESENTATIVE!

ASK ABOUT TRADE-IN CREDIT FOR YOUR OLD INVENTORY!

- “ The Prescott's OmniPlus/OSM200 is an excellent office microscope for the ENT surgeon. It is very affordably priced and functions exceptionally well. I give Prescott's my highest recommendation.”  
— Dr. Kuperan, *Houston Physicians Hospital*
- “ I like the ease of use. Everything is built in and I can record very easily.”  
— Yaara Berdan, *Herman Ostrow School of Dentistry of USC*
- “ Scope function is easy and simple. Visual acuity is top notch, and LED lighting is exactly what we needed.”  
— Dr. Montano, *McAllen Surgical Specialty Center*

Prescott's, Inc.  
18940 Microscope Way  
Monument, CO 80132

Phone: (800) 438-3937  
prescott@surgicalmicroscopes.com  
www.surgicalmicroscopes.com



EC REP  
EMERGO EUROPE  
Prinsessegracht 20  
2514 AP The Hague  
The Netherlands